




GV Health  
Life. Protected.



**Science-led  
Innovation**



**Smarter  
Solutions**

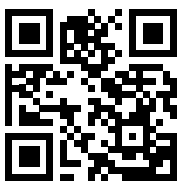

**Infection** Prevention

# Contents

<b>GV Health</b>	<b>4-14</b>
Our Mission	4
What we do	6
Where we work	6
Research & Innovation	8
Global Manufacture & Supply	10
Support & Education	12
GV App - Work. Smarter.	13
Sustainability	14
<b>Sharps Disposal</b>	<b>16-21</b>
<b>NEW</b> Life Protected Sharps Box	16-21
Sharps Box 5L (Yellow, Orange & Purple)	20
<b>Surfaces &amp; Environments - Cleaning &amp; Disinfection</b>	<b>22-33</b>
<b>NEW</b> Life Protected Cleaner	24-27
SoChlor Technology	28
SoChlor TAB - Disinfection - Multi-Purpose Disinfectant	30
SoChlor DST - Surface - Dual Clean & Disinfect	32
<b>Spill Kits - Containment, Removal &amp; Decontamination</b>	<b>34-43</b>
<b>NEW</b> Sorbent Technology & Spill Kits	36
<b>NEW</b> All-in-one Biohazard Spill Pack	38
Biohazard Spill Kits & Packs	38
Cytotoxic Spill Kit & Pack	40
Chemical Spill Kit & Pack	41
SoChlor Granules, High Absorbent Granules & UV Spill Kit	42-43
Life Protected Sorbent	43
<b>Polycare Healthcare Aprons</b>	<b>44-47</b>
Polycare Drug Round Tabards	48
Polycare Specimen Bags	50
<b>Polycare Healthcare Waste Bags</b>	<b>52-63</b>
Orange Clinical Bags	54
Yellow Clinical Bags	56
Tiger Stripe / Offensive Waste Bags	58
Cytotoxic Waste Bags	60
Black Waste Bags	62
Clear Waste Bags	62
Tags, Ties & Tape	64

Contract provider to the  
UK NHS for over 30 years



**GV Health**  
Life. Protected.






## Our Mission

# Creating safer environments for human health

Delivering integrated solutions, bringing together the best science and technology to the human challenges of Infection Prevention and environmental hygiene, through innovation and education.

**Life.** Protected.

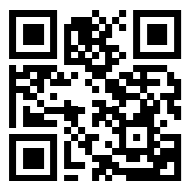



## What we do

In settings where infection prevention and environmental health are critical, our solutions save lives by achieving better patient outcomes in healthcare and improved public health.

We continuously invest in research and development to pioneer new product technologies and embrace more sustainable approaches. Integrated with our product offering, our education platforms and software systems, we provide powerful tools to help healthcare organisations further reduce their carbon footprint through data and information.

Our commitment to advancements in manufacturing and regulation have not only improved human health but helped reduce net impacts on the wider environment.


## Where we work

With a strong focus on infection prevention and control (IPC), we have established ourselves as a leader in the industry. Our waste segregation, environmental hygiene and safety solutions are trusted and used worldwide by healthcare professionals. From groundbreaking product design to comprehensive training and education, we offer unparalleled support to end-users in creating safer environments for human health.

We work closely with health providers around the world including the UK NHS where GV Health has been a proud contract provider for over 30 years.

## How we do this...

1. Our Research & Innovation
2. Our Global Manufacture
3. Our Knowledge and Support

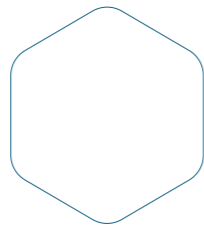




Life. Protected.

## Our Research & Innovation

From our brand new state-of-the-art technical facility in Stevenage, UK to our strategic partnerships around the world, we have developed and invested extensively in our resources and ability to deliver unique groundbreaking innovation.

As a Science and Technology first company, our unique capability to bridge scientific innovation with the hygiene needs of healthcare providers is pivotal in delivering impactful solutions. With a proven track record of leadership in formulating effective hygiene and infection prevention solutions, our research teams and state-of-the-art facilities, coupled with close collaborations with healthcare providers, allow us to fully understand their requirements and challenges. This synergy empowers us to pioneer industry-first solutions. While our global manufacturing network enables us to deliver these groundbreaking solutions at optimal value, addressing the budget challenges of today.

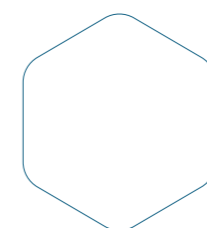

# Our Global Manufacture

To provide leading quality products, while meeting customer cost needs and demand, requires a comprehensive and global approach to manufacturing and development. Our international supply chain is central to delivering this, with our teams, partners and distribution channels connected around the world.

With manufacturing plants in the UK, Europe, the Middle East and Asia our robust and responsive supply chain has been pivotal in providing our solutions to create safe environments in which people live and work. While also ensuring continuity of supply through the increasing global challenges today.




# Our Knowledge Centre


## Insights. Learning. Support.

A core pillar of our solutions is rooted in our fundamental belief in the intrinsic importance of education and understanding for the end user. Providing industry leading products combined with our unique education and learning resources and platforms are critical in order to drive better outcomes in environmental hygiene. By providing end users and clinicians with a rich library of easy to follow content such as videos, articles and scientific insight, we can bring together and bridge the gap between end-user challenges and leading science, offering a more accessible insight and understanding of regulation and compliance.

*Improving the education around science and hygiene through dynamic impactful learning.*

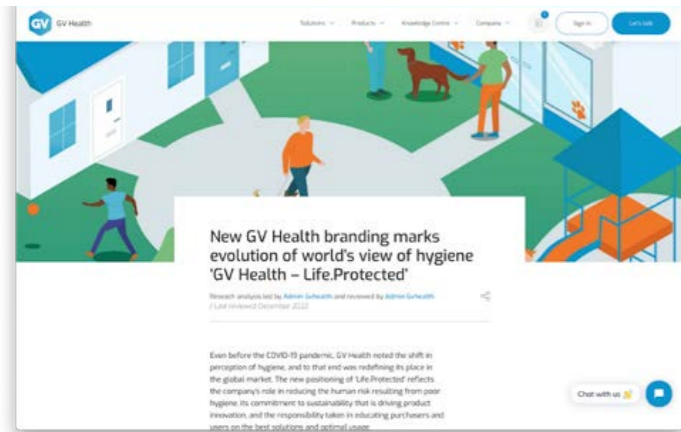

Insights & Articles

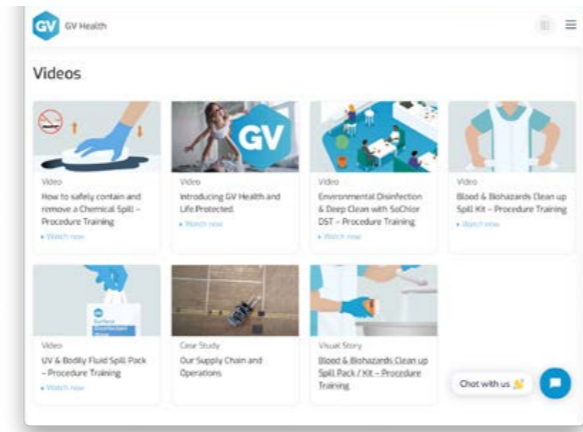

Videos

[gvhealth.com/knowledge-centre](http://gvhealth.com/knowledge-centre)




# GV App


## Reduce wastage. Reduce costs. Better Sustainability.

With the healthcare landscape rapidly evolving, organisations face two main responsibilities, environmental and economic. Maintaining low cost and carbon footprint can be challenging. The GV App offers comprehensive solutions for effective resource management, that balances environmental responsibility with financial efficiency.


Currently in use across the UK NHS in over 3,000 locations. Designed to help healthcare organisations track allocation of resources, reduce cost and improve sustainability by reducing wastage and improving efficiency. The GV App can help organisations save an average of 20% spend on their consumables. In addition, thanks to the smart data provided, it can also help monitor aspects of compliance across a healthcare organisation.

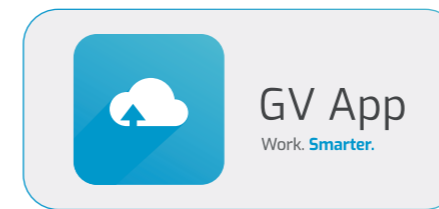

Developed in-house and in partnership with NHS Hospitals across the UK, the GV Health App can be used across all devices. Available via Browser, iOS and Android apps. The GV Health App is also used by a number of leading global Facility Management Companies or FMs.

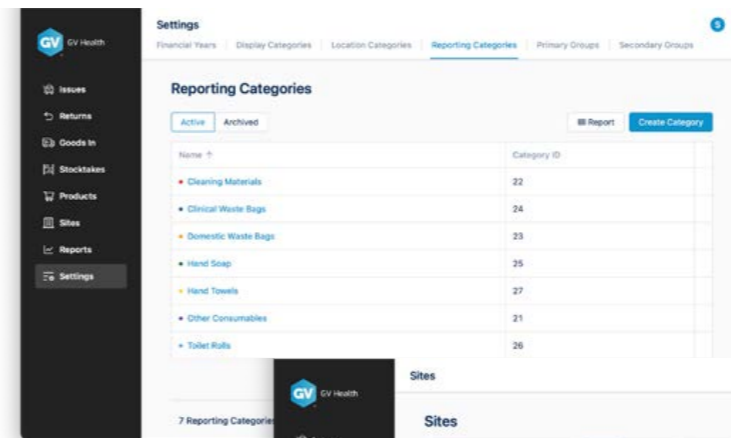

20% average saving per organisation

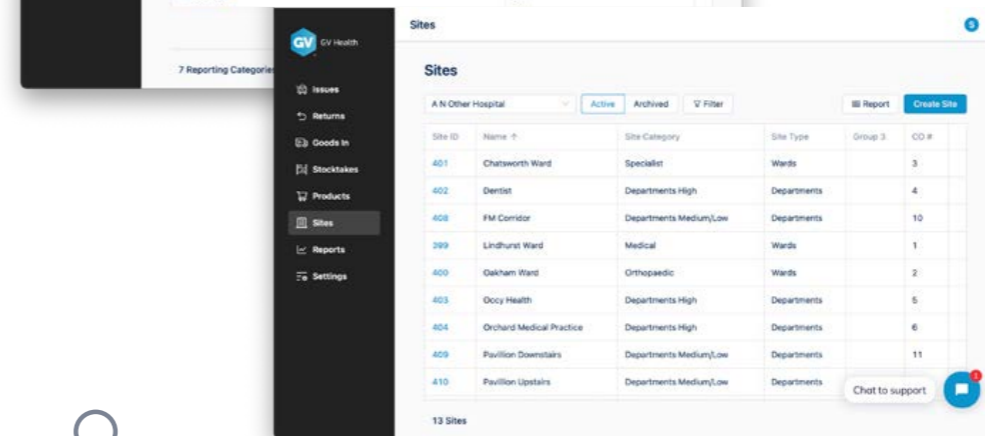

Available on iOS, Android & Web Browsers



### Request Trial or Demonstration

If you would like to learn more or request a demonstration, please contact your GV Health account manager.

[gvhealth.com/gv-app](http://gvhealth.com/gv-app)



## Sustainably Focused Since 1993

Protection is at the heart of GV Health. From our software to our technologies and innovations, we continue to drive more sustainable outcomes

### Our mission

Working alongside professionals committed to protecting human health and welfare, we see how vital it is that we also consider our impact on the planet and ensure we can provide cleaner, safer human environments now and in the future.

### Measurable and Transparent

Our commitment to sustainability continues to advance through all of our operations, optimising our existing product offerings, the reduction of virgin polymers in our bags, the introduction of new materials such as biopolymers to our products and the implementation of a Carbon Reduction Plan (CRP), Extended Producer Responsibility (EPR) outlining our accountability to sustainable practice, and our goals to achieve this. With a team of qualified experts including sustainability scientists, we are dedicated to exploring new technologies in infection control.

### Products, Technologies and Manufacture

Products that are more sustainable. Through technology, development, and design.

- Reduce manufacturing and delivery carbon footprint
- Reduce material resource consumption
- Reduce impact on environment at product end life



Our new Life Protected ranges of innovative and more sustainable products for the future of Infection Control and Prevention.

Life. Protected.

## 14001 accreditation



ISO 14001  
Certificate No. EM2001057

ISO 14001 demonstrates our Environmental Management Systems relating to the provision of our products and support services are maintained to the exacting standards required by ISO.

## Non-plastic Sharps Disposal

Developed with renewable biopolymers, this sustainable technology transforms the disposal of sharps.

Signifies a transition towards more sustainable choices in healthcare waste disposal

- Our cardboard flat-packed sharps boxes create a sustainable cycle by releasing biogenic carbon.
- Environmental impact of incineration is minimised in comparison to emissions associated with traditional fossil fuel-derived materials used in traditional plastic sharps boxes.
- Reduction in damaging emissions from incineration in comparison to plastic sharps boxes. Reducing the carbon footprint associated with incineration.

## GV App

Designed to track and monitor internal stock and material management, enabling real data analysis to identify and address inefficiencies to reduce wastage. The GV App is used in over 3,000 locations across the NHS. On average the system has helped reduced wastage by approx. 20%, substantially improving sustainability while also reducing costs.

The GV Health App, a software platform (pioneered for the NHS) for internal stock management, enabling real time monitoring and data analysis to identify and address inefficiencies to reduce material wastage (approx. 20% reduction in wastage achieved since its introduction in 1997), substantially improving sustainability.

**20%** reduction in usage of materials across multiple sites using the GV App

...meaning a **20%** improvement in sustainability of material use

## Life Protected Cleaner

82.4% reduction in carbon footprint of inner packaging compared to pre-concentrate cleaners.  
34% reduction in inner packaging waste (by weight) compared to pre-concentrate cleaners\*.  
Cardboard inner packaging reduces use of plastic.  
77% reduction in the carbon footprint and cost associated with transportation and storage compared to pre-concentrate cleaners  
Manufactured with 30% renewable energy.

## Reducing Plastic Since 1995

Optimised polymer blends to reduce waste bag thickness, thereby reducing plastic usage, since 1995, achieving up to a 76% reduction in material usage on core product lines. For our clinical bags, we have reduced the raw material used by 58 tonnes in our heavy-duty bags and 62 tonnes in medium-duty bags per year.

We are also working to increase the recycled content in our bags, thereby further reducing demand for virgin plastic, helping to further lower our carbon footprint.

**76%** reduction in material  
**120 tons** reduction per year in material used to make medium & heavy duty waste bags

**19%** reduction in packaging  
**30%** recycled polymer

## Sustainable IPC


SoChlor disinfectant tablet-based technology, designed to reduce transportation, packaging and storage requirements to reduce material wastage and carbon footprint.

The tablet delivery format of SoChlor reduces the need for transportation of water, creating lighter, more compact packages which create less CO2e during delivery. Further positives come from using SoChlor solutions with a reusable cloth, reducing the volume of disposable wipes and waste.



NEW

# NON-PLASTIC SHARPS BOX

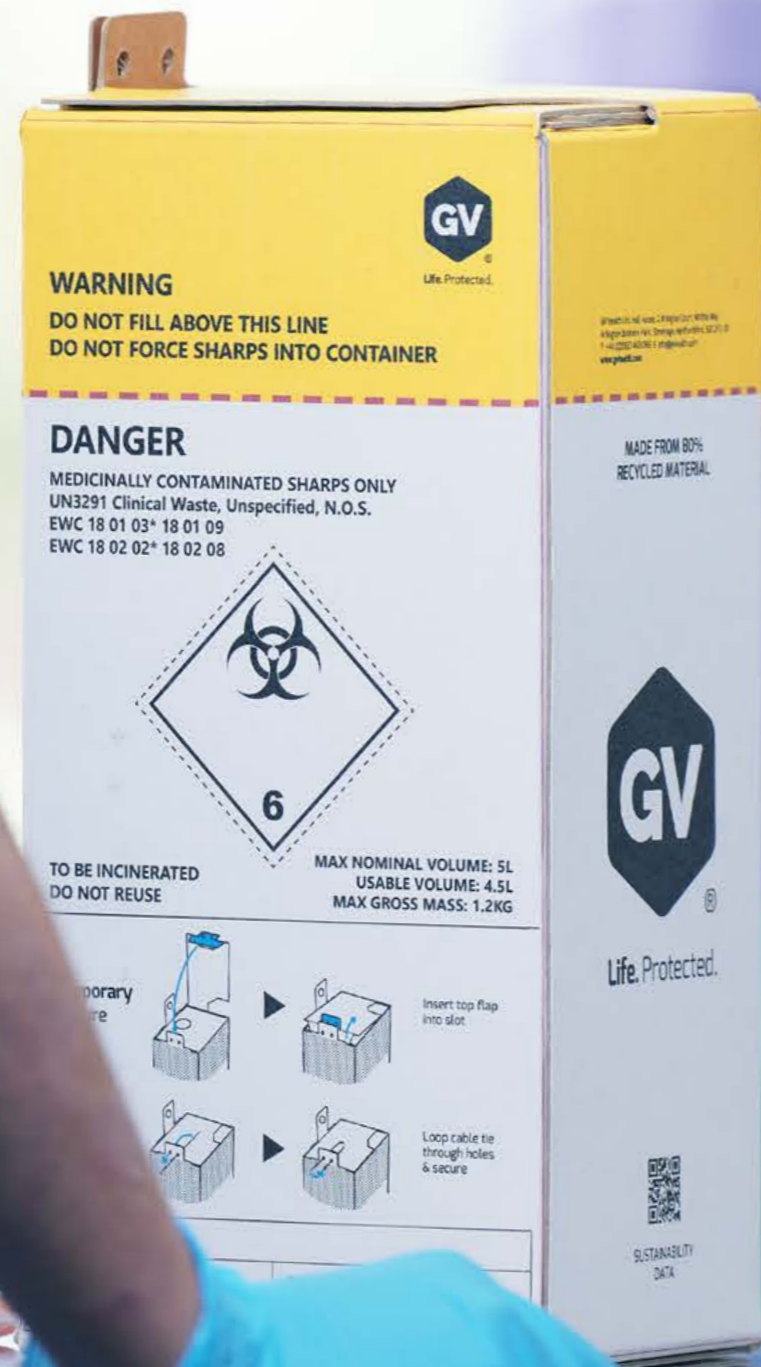

SHARPS BOX

SAFER  
SMARTER  
SUSTAINABLE



Life. Protected.

## Introducing our new Non-Plastic Sharps Box. The first of our Life Protected range of innovative & more sustainable IPC products

Developed with renewable biopolymers, this sustainable technology transforms the disposal of sharps.

GV Health's in-house custom developed cardboard composition ensures its unparalleled strength and integrity, making it a robust and reliable alternative to plastic.

Combined with our patent-pending two-part absorbent system comprised of a barrier and absorbent, our Sharps Box is capable of containing small amounts of liquids left over from syringes.

The flat-packed design significantly enhances transportation and storage efficiency. While its unique closure system ensures safe and secure containment.

Based on our research and available intel, we believe this is the first cardboard sharps box in the world to meet both BS EN ISO 23907-1 standards and ADR UN3291 certification.

# The future of sharps disposal...

## The first non-plastic Sharps Box available on the UK NHS framework

Our cardboard flat-packed sharps boxes are meticulously crafted through innovative research. Designed with a commitment to sustainability, its renewable cardboard and biopolymer-based liquid absorbent release biogenic carbon instead of fossil carbon when incinerated. Its flat-packed design is quick to assemble and optimises storage and shipping logistics fitting 1,350 boxes per pallet, unlike 432 of its plastic alternative, significantly lowering the carbon footprint. The cable-tie closure promotes secure disposal and joins together public health and sustainability, resonating with the future of healthcare. A forward-thinking solution, our sharps boxes are available in all waste streams and are a transformative leap in healthcare innovation.

- Cut down on incineration costs and environmental impact. Our cardboard flat-packed sharps boxes offer an environmentally-conscious alternative, decreasing the need for high-energy processes to incinerate traditional plastic sharps boxes.
- Production of biogenic carbon instead of fossil carbon.
- Flat-packed design for reducing carbon footprint reduction in transport.



### Easy assembly in 30 seconds

With a printed number system and easy-to-follow assembly instructions, healthcare professionals can set up the sharps boxes quickly and with ease. Enabling a hassle-free experience that combines speed, and user-friendly design, allowing for a seamless integration into existing waste management protocols.



### SHARPS BOX



Life. Protected.

#### PATENT PENDING

From design and technologies, to new manufacturing techniques and custom materials, the Life Protected Sharps Box has many first of their kind innovations

#### FULLY TESTED & COMPLIANT

✓ UN 3291 CERTIFIED PASSED ADR TEST REQUIREMENTS

✓ PASSED BS EN ISO 23907-1:2019 TEST REQUIREMENTS

- Tests including:
- Resistance to leakage
  - Resistance to needle penetration
  - Drop testing

### SAFER

#### Liquid leak prevention system

This unique feature, absent in other existing non-plastic sharps boxes, adds an extra layer of protection, setting our product apart as a comprehensive solution that prioritises not just safe disposal but also the containment of fluids. Crafted from biopolymer offers a sustainable alternative that reduces environmental impact.



#### Safety-first design

Our cardboard sharps boxes feature a cable-tied lock system designed to be permanent and secure, mitigating the risks of accidental openings and providing an extra layer of protection against sharp injuries, ensuring a secure healthcare environment.

Our purposeful lock system is engineered to be closed intentionally, eliminating the possibility of accidental closures, therefore reducing unnecessary waste of empty boxes. Embrace a safety-first approach without contributing to excess healthcare waste.



20% reduction transportation weight vs plastic

### SMARTER

#### Flat-pack design saves on shipping & storage

Flat-packed design enables efficient space utilisation, providing flexibility in storage without compromising on capacity. Store more units in less space in comparison to standard plastic sharps boxes, allowing allocation of more space for other equipment, streamlining inventory management within hospitals and clinics.

Optimised shipping logistics: By accommodating more boxes per pallet compared to standard plastic sharps boxes, the Life Protected Sharps Box significantly reduces the carbon footprint associated with transportation.

When compared to existing 5 Litre Plastic Sharps Boxes on the market today:


Existing Plastic Sharps Box: 432 per pallet with a weight of 365g per box  
Life Protected Sharps Box: 1,350 per pallet with a weight of 300g per box

85% less space & Storage vs Plastic Boxes



### SUSTAINABLE

#### Biogenic Closed-Loop Carbon System

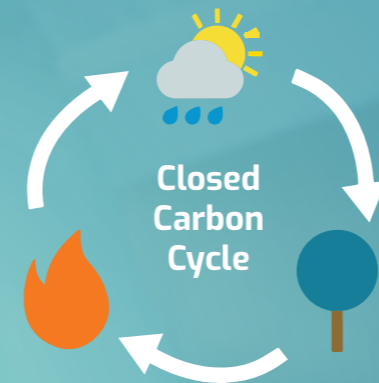

Signifies a transition towards more sustainable choices in healthcare waste disposal


- Our cardboard flat-packed sharps boxes create a sustainable cycle by releasing biogenic carbon.
- Environmental impact of incineration is minimised in comparison to emissions associated with traditional fossil fuel-derived materials used in traditional plastic sharps boxes.
- Reduction in damaging emissions from incineration in comparison to plastic sharps boxes. Reducing the carbon footprint associated with incineration.

✓ Made from FSC certified cardboard

First Non-plastic Sharps Box on the NHS Supply Chain

#### Smart Wall Mount System

The Smart Wall Mount System is fast and easy to secure with a custom adhesive that holds in place on any surface. Can also be easily removed via release mechanism without leaving marks.

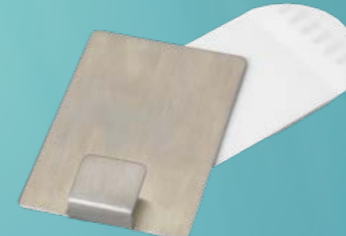

# Sharps Box - Yellow

**Code: LPSCY5** **ISO 23907-1 Passed** **ADR Certified**

For the disposal of partially discharged sharps including those contaminated with medicines other than those that are cytotoxic and cytostatic.

Comes as an easy-to-assemble flat pack.

**PACK CONTENTS:**  
15 x Yellow Sharps Box  
15 x Zip-Tie

**UN 3291 CERTIFIED**  
**PASSED ADR TEST REQUIREMENTS**

**PASSED BS EN ISO 23907-1:2019**  
**TEST REQUIREMENTS**

Tests including:


- Resistance to leakage
- Resistance to needle penetration
- Drop testing

✓ **Made from**  
FSC certified cardboard

**Clinical Waste Disposal**

MEDICINALLY CONTAMINATED SHARPS ONLY  
UN3291 Clinical Waste, Unspecified, N.O.S.  
UN 4G/Y2/S/24/GB/7568

EWC 18 01 03\* and 18 01 09  
EWC 18 02 02\* and 18 02 08

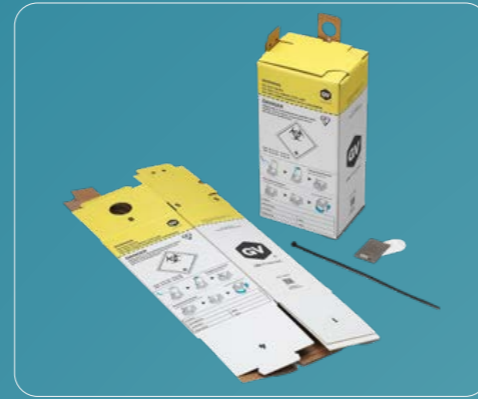

Pack Code	Product	Size	Quantity	Availability
LPSCY5	Life Protected Sharps Box - Yellow	5 litre	15 per pack	Q1 2025

# Sharps Box - Purple

**Code: LPSCP5** **ISO 23907-1 Passed** **ADR Certified**

For the disposal of sharps including those contaminated with cytotoxic and cytostatic medicines.

Comes as an easy-to-assemble flat pack.

**PACK CONTENTS:**  
15 x Purple Sharps Box  
15 x Zip-Tie

**UN 3291 CERTIFIED**  
**PASSED ADR TEST REQUIREMENTS**

**PASSED BS EN ISO 23907-1:2019**  
**TEST REQUIREMENTS**

Tests including:

- Resistance to leakage
- Resistance to needle penetration
- Drop testing

✓ **Made from**  
FSC certified cardboard

**Clinical Waste Disposal**

MEDICINALLY CONTAMINATED SHARPS ONLY  
UN3291 Clinical Waste, Unspecified, N.O.S.  
UN 4G/Y2/S/24/GB/7568

EWC 18 01 03\* and 18 01 08\*  
EWC 18 02 02\* and 18 02 07\*



Pack Code	Product	Size	Quantity	Availability
LPSCP5	Life Protected Sharps Box - Purple	5 litre	15 per pack	2025

# Sharps Box - Orange

**Code: LPSC05** **ISO 23907-1 Passed** **ADR Certified**

For the disposal of sharps that are not contaminated with medicinal products.

Comes as an easy-to-assemble flat pack.

**PACK CONTENTS:**  
15 x Orange Sharps Box  
15 x Zip-Tie

**UN 3291 CERTIFIED**  
**PASSED ADR TEST REQUIREMENTS**

**PASSED BS EN ISO 23907-1:2019**  
**TEST REQUIREMENTS**

Tests including:

- Resistance to leakage
- Resistance to needle penetration
- Drop testing

✓ **Made from**  
FSC certified cardboard

**Clinical Waste Disposal**

MEDICINALLY CONTAMINATED SHARPS ONLY  
UN3291 Clinical Waste, Unspecified, N.O.S.  
UN 4G/Y2/S/24/GB/7568

EWC 18 01 03\* and 18 02 02\*



Pack Code	Product	Size	Quantity	Availability
LPSC05	Life Protected Sharps Box - Orange	5 litre	15 per pack	2025



## REGISTER TO ENROL FOR FIRST STOCK AVAILABILITY

First stock will be available Q1 2025. Register interest with your Account Manager or contact us directly



# ENSURING SAFER ENVIRONMENTS



## Smart, advanced and forward looking technologies for cleaning and disinfection

Our surface cleaning and disinfectant solutions have been developed to meet the most difficult challenges in healthcare infection prevention.

Designed to provide the best outcomes on floors and hard surfaces.

Our products are designed to flex and adapt to the current needs of the environment.

Developed in healthcare and now used across sectors to provide the highest standards of hygiene.





**Life Protected Cleaner**  
Water Soluble Matrix

**2024**  
**Innovative Product**  
**of the Year**  
Awarded  
**Highly Commended**  
**Cleaning Excellence Awards**



SCAN TO LEARN  
MORE ABOUT  
LPC


**Innovative Product of the Year 2024**  
 Awarded Highly Commended  
 Cleaning Excellence Awards

**Life Protected Cleaner**  
 Water Soluble Matrix

**The future of cleaning technology**

Introducing a breakthrough in cleaning technology with GV Health's new Life Protected Cleaner (LPC). LPC water-soluble sheets, dissolve in water in seconds to create a high-performance cleaning solution for environmental surfaces. Developed with an innovative surfactant formulation based on GV Health's new Water Soluble Matrix Technology (WSM) which is set to revolutionise the cleaning industry.

**Simply dissolve sheets in water...**

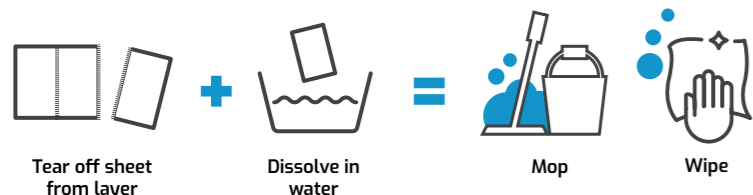

36 sheets inside each pack



**Smarter**

**Easy dosing**

With the cleaning formula already pre-measured in each LPC sheet, users no longer have to worry about measuring and pouring the right amount. This saves time and eliminates the risk of over-use or under-use, which can result in ineffective cleaning or waste. **Dissolves at least 10x faster** than all other solid surface cleaning products, even in cold water. Safe and easy to dispense – no splash risk with the solid sheet and no need for manual pump dispensers. Suitable for use with cleaning machines.

**High performance**

This multi-purpose cleaning formulation enables a high level of cleaning, removing dirt, grease and tacky residues, saving both time and labour costs.

**At least x4 less storage space**

LPC sheets are easy to use and store, making them a convenient choice for people on-the-go. LPC sheets reduce storage space required by at least 4 times vs liquid concentrate cleaners.

**Versatile use & application**

Compatible for use with spray bottles, spray, jet or foaming nozzles, automated cleaning equipment, and mops, the LPC sheets can be used with a broad range of application methods.

**Multiple applications**

LPC's multiple applications mean that it can easily replace a number of existing cleaning products, making it a more storage efficient, cost effective and sustainable alternative.



**Up to 90L of solution from 1 pack**

LPC's compact design means that one pack can make up to **36L** of ready-to-use solution for hard surface cleaning and **90L** for floor cleaning.

Application	Usage rate		Total ready to use solution volume per SKU (litres)	
	Litres per sheet	Sheets per litre	Pack	Case
<b>Hard Surfaces</b>	1	1	36 L	<b>1,728 L</b>
<b>Floors</b>	2.5	0.4	90 L	<b>4,320 L</b>

**The Science & Technology**

LPC is a revolutionary cleaning technology, developed in GV Health's state of the art laboratory facilities. Based on new chemistries, test methods and manufacturing processes, LPC is designed to meet the complex, varied and rapidly evolving needs of facility hygiene operations.

**Water Soluble Matrix**

LPC is based on GV Health's Water Soluble Matrix (WSM) technology platform.

- Up to **50%** reduction in chemical use compared to other solid cleaning products.
- Reduced PPE - Safer formula means only gloves are required
- Safer and easier to handle due to a **95.5%** reduction in weight compared to liquid pre-concentrate cleaners.
- A safer, pH neutral surfactant-based formula, kinder to materials and compatible with a wide range of surfaces.
- Leaves no surface residues, unlike most solid effervescent cleaning products.

**Significantly more Sustainable**

- **82.4%** reduction in carbon footprint of inner packaging compared to pre-concentrate cleaners.
- **34%** reduction in inner packaging waste (by weight) compared to pre-concentrate cleaners\*.
- Cardboard inner packaging reduces use of plastic.
- **77%** reduction in the carbon footprint and cost associated with transportation and storage compared to pre-concentrate cleaners.
- Manufactured with **30% renewable energy**.

Code	Product	Pack Size	Dilution / Use
<b>LPC036</b>	Life Protected Cleaner	Pack of 36 sheets (12 layers, 3 sheets per layer)	1 sheet in 1 litre for hard surfaces 1 sheet in 2.5 litres for floors
<b>LPCCSB</b>	Life Protected Spray Bottle	Unit	750ml



**Hard Surfaces**  
 Wipe (Foam, Spray or Soak)



**Floors**  
 Mop or Machine compatible

Developed for use in healthcare and professional cleaning, LPC is already used in the NHS, HMP, hospitality and other sectors to ensure the best outcomes in hygiene and cleaning.



# Unique Technology

## Why SoChlor is one of the most advanced antimicrobials on the market today

Used around the world, SoChlor is one of the most effective and versatile disinfectant and decontamination products on the market today. It plays a central role in protecting human health everyday across sectors, by ensuring best practice hygiene and breaking chain of infection from harmful organisms.

From high level disinfection for infection control and prevention in healthcare, to a broad range of hygiene applications in agriculture, catering, hospitality, commercial and home use; SoChlor is trusted by professionals and users around the world, to achieve the highest standards

and ensure the safest outcomes for human health wherever it is used.


SoChlor is a tablet-based solution offering high efficacy disinfection and decontamination, while being highly versatile. At the heart of SoChlor is a Chlorine based chemistry – hypochlorous acid which provides industry leading performance, with up to 100 times the anti microbial effectiveness of other chlorine products such as bleach, while providing a range of further benefits and advantages over many actives in the market.

## The technology and benefits of SoChlor

SoChlor achieves higher performance and greater benefits over other products, to deliver proven, reliable, and trusted performance.

### Key benefits include:

1. High Antimicrobial performance
2. Smart Tablet Format
3. Intelligent Active Delivery
4. Inspired by Nature
5. Versatile applications across multiple sectors

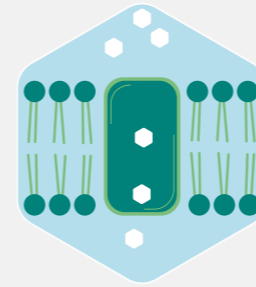


## High Antimicrobial Performance

High Antimicrobial performance. SoChlor technology is effective at killing a wide range of pathogens, including Bacteria, Viruses, Spores, and Fungi. Particularly pathogens of high concern such as, mycobacteria, non-enveloped viruses and spore forming bacteria. And can be used for high-level disinfection in hospitals, helping to break the chain of infection, interceding at the early stages of transmission. SoChlor's antimicrobial effectiveness is thanks to the ability of Hypochlorous acid, which unlike other antimicrobials, is able to pass through the pathogen's cell membrane, thereby entering the cell. Cells are designed to protect themselves from the external environment, but not from within. Therefore, as an oxidising agent, once within the cell, Hypochlorous acid is able to more effectively cause cell death. This enables Hypochlorous acid to act as a high-level disinfectant when used as directed.

## Smart Delivery

SoChlor technology provides smart delivery in two key ways:


**Tablet:** Tablet format doesn't contain water, minimising the need for packaging, transportation and storage space, thus reducing carbon footprint compared to water based products.

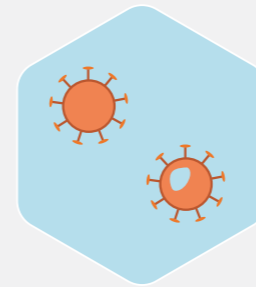

**Active:** When the NaDCC in the SoChlor detects the presence of pathogens it releases more Hypochlorous acid. How does it achieve this? Thanks to the reversible reaction of NaDCC with water, it creates an equilibrium, which means that firstly: only a portion of the total available Hypochlorous acid is delivered at any given time and secondly, more Hypochlorous acid is delivered if it comes in to contact with pathogens.

## Safe and naturally occurring active


SoChlor is biomimetic, taking inspiration from nature in action, developed from millions of years of evolution. This is because the Hypochlorous acid produced by NaDCC is similar in some respects to the immune response of the Human body, which also releases Hypochlorous acid at the site of infection to combat invading pathogens. Given its prevalence in nature Hypochlorous acid has been described as Nature's own powerful non-toxic disinfectant.

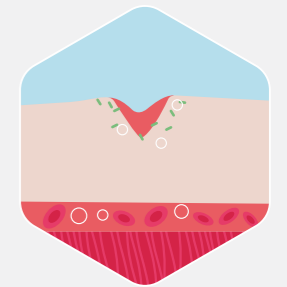

## Versatile use and applications across sectors

SoChlor is one of the most versatile antimicrobial products for disinfection on the market today, to protect human health.

- Healthcare: general disinfection, outbreak disinfection/decontamination, biohazard spill disinfection / decontamination, daily clean, deep clean, terminal clean and isolation clean.
- Decontamination in clinical and laboratory settings: stainless steel instruments, and Laboratory discard jars.
- Disinfection of baby and infant equipment: infant bottles and Teats.
- Catering and food: disinfection of food preparation surfaces, vegetable and salad washing.
- Domestic and commercial cleaning: WC's, drains, sinks, dish cloth and mops.
- Agriculture: disinfection of drinking water, disinfection of animal housing floors and pens, bedding, surfaces and walkways, equipment, footwear baths, and vehicles.

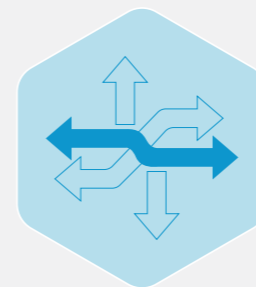

# SoChlor<sup>®</sup> TAB

## Multi-Purpose Disinfectant Tablets

An essential in infection prevention and control (IPC) across sectors, SoChlor (Sodium dichloroisocyanurate, NaDCC) is a smart, sustainable precursor to hypochlorous acid – a chlorine-based, safe and powerful solution to disinfection.

## High antimicrobial performance

With a portfolio of ready-to-use tablets, SoChlor TAB tablets dissolve to formulate the correct concentration of disinfectant, every time - simplifying disinfection processes and eliminating product waste.

SoChlor's biomimetic technology emulates the same response as certain white blood cells in the body, where a non-toxic, biodegradable compound hypochlorous acid is released on exposure to pathogens. Activated upon dissolving in water, NaDCC releases hypochlorous acid – the safe, smart active in SoChlor. Holding a much higher level of available chlorine, hypochlorous acid is approximately 100x more effective than household bleach, meaning SoChlor's anti-microbial performance is efficacious against more resistant pathogens including mycobacteria, spore-forming bacteria and non-enveloped viruses. SoChlor is therefore suitable for high-level and outbreak disinfection, such as against SARS-CoV-2.

SoChlor TAB's stable, tablet formula and small composition makes for easy transportation. With higher quantities per delivery due to reduced packaging, they are a more sustainable option for infection control.

## Tested & certified:

SoChlor<sup>™</sup> TAB tablets are tested and certified against EN test standards: B5:EN1650; B5:EN13704; B5:EN1276 (EN13727); B5:EN14476 (Adenovirus); B5:EN14476 (Poliovirus); B5:EN14476 (Norovirus); B5:EN14348; B5:EN13624; B5:EN13727; B5:EN13697 and conform to UK Department of Health, HPA and WHO guidelines for isolation, terminal and general environmental cleaning and disinfection.



SCAN TO WATCH SOCHLOR TECHNOLOGY VIDEO

SOCHLOR TAB



## Versatile use and applications across sectors

### Healthcare

SoChlor TAB tablets are a versatile disinfecting solution across all healthcare settings. For general, multi-area disinfection, SoChlor TAB can be used on surfaces across wards, labs, clinics and in theatre, diluted at 1,000 ppm as recommended by the Department of Health and the Health Protection Agency. It is also effective for the disinfection of laboratory instruments including laboratory discard jars and stainless steel instruments, and at higher concentrations, can be used for outbreak disinfection. At 10,000 ppm, SoChlor TAB tablets are effective for the disinfection of biohazardous spills and bodily fluid contamination.



### Primary and Social Care

Highly effective surface disinfection makes SoChlor an ideal for use in care facilities, particularly on high contact equipment including bed frames, commodes and domestic equipment. Use in communal areas, kitchens and bathrooms is ideal to effectively disinfect surfaces contaminated with bodily fluids, including urine and vomit, killing pathogens to subsequently reduce risk of infection in high-risk patient settings.



### Agriculture & Veterinary

Suitable for the high-level disinfection of instruments and the disinfection of biohazardous materials, SoChlor's profile of versatile applications make it an optimal solution for the sanitation and disinfection of clinical veterinary settings, as well as in general areas such as waiting rooms and offices. SoChlor technology kills pathogens and reduces risk of cross-infections, providing clean and safe environments for animals in your care. Agricultural applications of SoChlor<sup>™</sup> include the disinfection of feeding equipment, shelters and animal housings, and is safe enough to be used for the purification and sanitation of water.



### Education & Childcare

SoChlor's technology means it is suitable for general and surface disinfection across a range of settings, including across classrooms, halls and communal areas. When used according to instructions, SoChlor is safe enough to be used in paediatric health and social care settings, to disinfect equipment, toys and even infant bottles and teats.



## ORDER INFORMATION

A wide range of tablet sizes to best fit your use.

Order Code	Product	NaDCC / Chlorine per tablet	Quantity
<b>MFB694</b>	SoChlor TAB 4.5 100	4.5g / 2.5g	100 Tablets per pot / 6x pots per case
<b>MFB255</b>	SoChlor TAB 2.5 100	2.5g / 1.5g	100 Tablets per pot / 6x pots per case
<b>MFB637</b>	SoChlor TAB 2.5 200	2.5g / 1.5g	200 Tablets per pot / 6x pots per case
<b>MFB254</b>	SoChlor TAB 1.7 200	1.7g / 1.02g	200 Tablets per pot / 6x pots per case
<b>MFB605</b>	SoChlor TAB 0.5 600	0.5g / 0.3g	600 Tablets per pot / 6x pots per case



# SoChlor<sup>®</sup> DST

## Disinfectant, Surfactant Tablets

SoChlor DST Disinfectant Tablets clean and disinfect in one process. The combined action of the surfactant, chlorine and oxidising action deeply cleans and decontaminates.

### Dual-action: Clean & Disinfect

The dual action Disinfectant Surfactant Tablets or "DST" deliver effective environmental cleaning and surface disinfection for healthcare and public health.

SoChlor (Sodium dichloroisocyanurate, NaDCC) is a smart, sustainable precursor to hypochlorous acid – a chlorine-based, safe and powerful solution to disinfection. Combined with a compatible detergent – have the ability to clear and disinfect in one simple step. SoChlor DST works to remove light surface loads, followed with disinfecting action, to sanitise and disinfect surfaces quickly. SoChlor DST tablets dissolve to formulate the correct concentration of solution, every time - eliminating guess work and reducing user error – resulting in hassle-free cleaning and disinfection while eliminating product waste.

SoChlor's biomimetic technology emulates the same response as certain white blood cells in the body, where a non-toxic, biodegradable compound hypochlorous acid is released on exposure to pathogens. Activated upon dissolving in water, NaDCC also releases hypochlorous acid – the safe, smart active in SoChlor. Holding a much higher level of available chlorine, hypochlorous acid is approximately 100x more effective than household bleach, meaning SoChlor's anti-microbial performance is efficacious against more resistant pathogens including mycobacteria, spore-forming bacteria and non-enveloped viruses. SoChlor™ is therefore suitable for high-level and outbreak disinfection, such as against SARS-CoV-2.

SoChlor DST pairs SoChlor's high anti-microbial action with a compatible detergent, DST transforms the cleaning, rinsing and disinfecting process into one quick step. Its stable, tablet formula and small composition makes for easy transportation. With higher quantities per delivery due to reduced packaging, SoChlor DST tablets are a more sustainable option for cleaning and infection control.



SCAN TO WATCH SOCHLOR TECHNOLOGY VIDEO

SOCHLOR DST



### Tested & certified

SoChlor DST tablets conform to UK Department of Health, HPA and WHO guidelines for isolation, terminal and general environmental cleaning and disinfection. They are also tested and certified against EN test standards:

**EN 1276, EN 13697, EN 13727, EN 16615, EN 17126, EN 13624, EN 14348, EN 16615, EN 14476**

**Effective against Mycobacteriacidal, Bacterial, Fungal, Yeastocidal, Virucidal & Sporocidal Activity**

### Keeping surfaces safe during COVID-19

Chosen as a key product for the UK Government's COVID-19 response, SoChlor DST represented more than 50% of the UK's environmental hygiene solutions in settings including NHS hospitals, care homes, prisons, testing & vaccination centres and other key government facilities.

### Fast. Effective. Trusted. With leading support

When it comes to Infection Prevention and Control (IPC), we believe the product is only half the story to create effective solutions. Support and training are essential to ensure best practice and implementation to achieve the highest hygiene standards for IPC.



### High-level Disinfection


- Safe for daily disinfection of wards and hospital settings
- Efficacious against bodily fluids, including blood, making it effective against spills
- Powerful for use during outbreaks of infection

### Applications

- Terminal, isolation and outbreak cleaning.
- Environmental cleaning of rooms, bays and wards
- Cleaning of bed spaces and patient area surfaces
- Cleaning of toilet and bathroom areas
- Mattresses covers
- Enamelled surfaces
- Stainless steel
- Sealed floors
- Worktops and other hard surfaces



Scan to watch the DST video guide

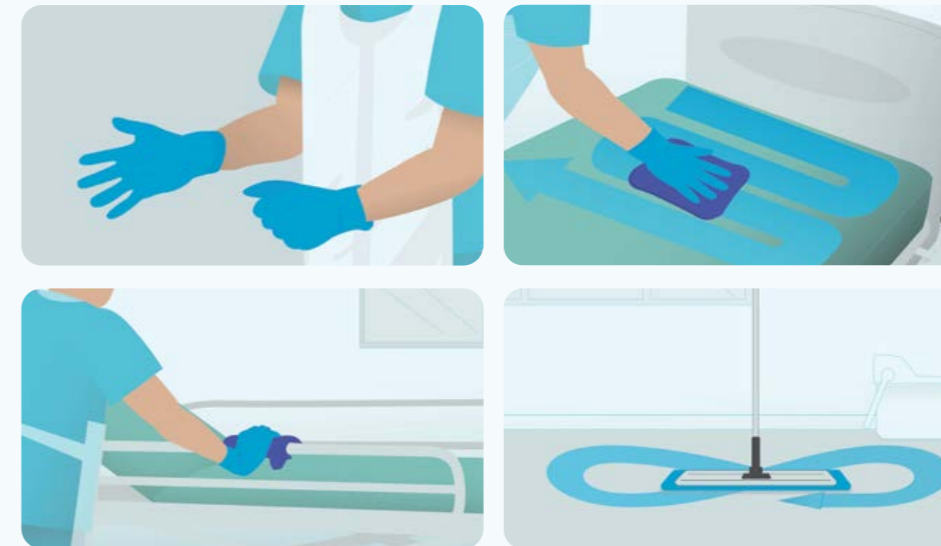

### ORDER INFORMATION

Order Code	Product	NaDCC / Chlorine per tablet	Quantity
<b>MFB256</b>	SoChlor DST 1.7 100 (Dual Clean & Disinfect)	1.7g / 1.02g	100 Tablets per pot / 6x pots per case
<b>MFB294</b>	SoChlor DST 1.7 200 (Dual Clean & Disinfect)	1.7g / 1.02g	200 Tablets per pot / 6x pots per case
<b>MFB257</b>	SoChlor Diluter 1 Litre		
<b>MFB258</b>	SoChlor Diluter 2 Litre		
<b>MFB259</b>	SoChlor Diluter 5 Litre		
<b>MFB693</b>	SoChlor Spray Diffuser 1 Litre		

# SPILL KITS & PACKS



## Health & Safety Spill Management for all sectors

Used across all sectors to ensure the highest standards of IPC and health and safety.

GV Health spill kits meet the highest efficacy and personal protection standards specified by the World Health Organisation, the UK Health and Safety Executive, UK Department of Health and CDC.



Scan to learn more about our biohazard spill kits



# Introducing our new Life Protected sorbent technology for our Spill Kits & Packs

Discover the innovation of GV Health's LP Multipurpose Sorbent Powder - a revolutionary spill management solution designed for compatibility with follow-on chlorine-based disinfectants, ensuring the surface is safely and effectively disinfected afterwards. The LP Sorbent Powder efficiently absorb and adsorb water, polar non-aqueous, and non-polar liquids. With superior containment properties and a hazard-free formula, they enable safe usage on a range of surfaces. Experience high-performance sorption, coupled with acid and urine neutralisation, heavy metal binding, and deodorising action, streamlining clean-up for all types of spills. The intelligent, non-hazardous, and sustainable design ensures rapid, safe disposal, making it an important tool for effective infection control.

### Multipurpose Sorbent Powder

A first-to-market and compatible to use with chlorine disinfecting chemicals, our granules redefine the landscape for spill management. With a sorbent technology optimised to both contain and neutralise a vast range of fluids, experience a new standard for safer and more effective surface disinfection.

### Unrivalled Versatility

High-Performance sorption: Elevate your spill response with our powder, unmatched in retaining more liquid on a weight/weight basis than any other mineral absorbent on the market.

Unlike other granules on the market, GV Health's Multipurpose Sorbent Powder is able to absorb:

- Water
- Polar non-aqueous liquids

**Non-polar liquid:** Maintains a consistent sorption rate for more efficient containment. Effectively absorbs and adsorbs water-based solutions without precipitating their dissolved contents, keeping spills properly contained.

**Acid & Urine Neutralisation:** One-product, multi-spill solution - chemically neutralises acidic solutions (up to 0.1M monoprotic acidic solutions), including human vomit and urine (by binding ammonium cations).

**Heavy Metal Binding:** Capable of chemically binding heavy metal cations, including Barium, Lead, Nickel, Arsenic, Mercury, Cadmium, Caesium, Copper, and Chromium.

**Deodorising action:** Experience a new standard of cleanliness with our powder, which not only physically absorb and adsorb but also chemically bind volatile organic molecules, eliminating unpleasant odours.

### Zero hazards & enhanced safety:

- Formulated without aggressive chemicals, so suitable for use across a wide range of surfaces, including carpets and fabrics.
- Minimal expansion during sorption, allowing spills to be properly contained, and reducing choking hazards if accidentally ingested.
- Non-hazardous formulation and minimal dusting during use.
- Creates a non-slip surface.
- No toxic gas generation during incineration.
- Safe for use on acidic spills.
- Clumps together upon sorption of water based liquids, making it easier to collect with (cardboard) scoop and scraper and minimises cross-contamination and risk of healthcare-associated infections (HCAIs).

**Environmentally Sustainable:** Contains natural minerals and polymers and is free from petroleum-based products.

**Longer Shelf Life:** Specially developed manufacturing process for greater stability, ensuring components remain homogenous during storage.



Life Protected Sorbent Powder

**NEW**

## ALL-IN-ONE Biohazard Spill Pack FOR ALL BIOHAZARD SPILLS

**LP Sorbent Powder**

- MJZ005 Pack contents:**
- 1 x pair of nitrile gloves
  - 1 x apron
  - 1 x **NEW LP Sorbent Powder Sachet**
  - 1 x diluter
  - 1 x **NEW cardboard scoop and scraper**
  - 1 x clinical waste bag
  - 4 x SoChlor TAB tablets
  - 1 x **NEW SoChlor Applicator**



**PLUS UPDATES TO OUR ENTIRE RANGE**  
 #NEW DESIGNS  
 #NEW NEXT GENERATION PERFORMANCE  
 #MORE SUSTAINABLE



Scan to to learn more about LP Sorbent

## Biohazard Spill Pack (All-in-one)

**Code: MJZ005** **LP Sorbent Powder** **All Biohazard Spills**

Enabled by the latest sorbent innovation and technology, the new MJZ005 Biohazard Spill Pack all-in-one can be used on all Biohazard Spills while also providing the 10,000ppm recommended by the WHO and healthcare organisations around the world for spills such as blood.

**PACK CONTENTS:**

- 1 x pair of nitrile gloves
- 1 x apron
- 1 x New LP Sorbent Powder Sachet
- 1 x diluter
- 1 x New cardboard scoop and scraper
- 1 x clinical waste bag
- 4 x SoChlor TAB tablets
- 1 x New SoChlor Applicator

**SPILLS COVERED INCLUDE:**

- |                           |                    |
|---------------------------|--------------------|
| <b>Blood</b>              | Semen              |
| Amniotic fluid            | Breast Milk        |
| Unfixed tissues/organs    | Vaginal secretions |
| Peritoneal fluid          | Saliva             |
| Pleural fluid             | <b>Urine</b>       |
| Pericardial fluid         | <b>Vomit</b>       |
| Synovial fluid            | Faeces             |
| Cerebrospinal fluid (CSF) | Sputum             |

**NEW**



## Biohazard Spill Pack

**Code: MJZ003** **SoChlor Granules**

**PACK CONTENTS:**

- 1 x pair of nitrile gloves
- 1 x apron
- 1 x 100g SoChlor NaDCC Granules
- 1 x diluter
- 1 x New cardboard scoop and scraper
- 1 x clinical waste bag
- 4 x SoChlor TAB tablets
- 1 x New SoChlor Applicator

**SPILLS COVERED INCLUDE:**

- |                        |                           |
|------------------------|---------------------------|
| Blood                  | Cerebrospinal fluid (CSF) |
| Amniotic fluid         | Semen                     |
| Unfixed tissues/organs | Breast Milk               |
| Peritoneal fluid       | Vaginal secretions        |
| Pleural fluid          | Saliva                    |
| Pericardial fluid      |                           |
| Synovial fluid         |                           |



## UV Spill Pack

**Code: MJZ004** **High Absorbency Granules**

**PACK CONTENTS:**

- 1 x pair of nitrile gloves
- 1 x apron
- 1 x sachet of High Absorbency Granules
- 1 x New cardboard scoop and scraper
- 2 x packs of disinfectant surface wipes
- 1 x clinical waste bag

**SPILLS COVERED INCLUDE:**

- Urine
- Vomit
- Faeces
- Saliva
- Sputum



## Biohazard Spill Kit Mini - 2 pack

**Code: MJZ016** **Refill Pack Format** **All Biohazard Spills**

Designed to cover all types of Biohazard and Body Fluid Spills. Small compact design allows portability and storage anywhere to ensure easy accessibility.



**BIOHAZARD SPILL PACK X1**

- Refill pack code: MJZ003
- 1 x pair of nitrile gloves
  - 1 x apron
  - 1 x 100g SoChlor NaDCC Granules
  - 1 x diluter
  - 1 x New cardboard scoop and scraper
  - 1 x clinical waste bag
  - 4 x SoChlor TAB tablets
  - 1 x New SoChlor Applicator

**URINE/VOMIT SPILL PACK X1**

- Refill pack code: MJZ004
- 1 x pair of nitrile gloves
  - 1 x apron
  - 1 x sachet of High Absorbency Granules
  - 1 x New cardboard scoop and scraper
  - 2x sachets of IPA surface wipes
  - 1x clinical waste bag
  - Illustrated instructions on pack



## Biohazard Spill Kit Duo - 2 pack

**Code: MJZ010** **Refill Pack Format** **All Biohazard Spills** **LP Sorbent Powder**

Designed to cover all types of Biohazard and Body Fluid Spills. Small compact design allows portability and storage anywhere to ensure easy accessibility.

**NEW**

**BIOHAZARD SPILL PACK(ALL-IN-ONE) X2**

- Refill pack code: MJZ005
- 1 x pair of nitrile gloves
  - 1 x apron
  - 1 x New LP Sorbent Powder Sachet
  - 1 x diluter
  - 1 x New cardboard scoop and scraper
  - 1 x clinical waste bag
  - 4 x SoChlor TAB tablets
  - 1 x New SoChlor Applicator

**SPILLS COVERED INCLUDE:**

- |                           |                    |
|---------------------------|--------------------|
| <b>Blood</b>              | Semen              |
| Amniotic fluid            | Breast Milk        |
| Unfixed tissues/organs    | Vaginal secretions |
| Peritoneal fluid          | Saliva             |
| Pleural fluid             | <b>Urine</b>       |
| Pericardial fluid         | <b>Vomit</b>       |
| Synovial fluid            | Faeces             |
| Cerebrospinal fluid (CSF) | Sputum             |



Available Q1 2025

## Biohazard Spill Kit Midi

**Code: MJZ018** **SoChlor Granules**

For use on up to 10 Biohazard Spills.



**MIDI 10 SPILL KIT CONTENTS:**

- 100 x 2.5g SoChlor TAB tablets
- 1 x 500g SoChlor NaDCC Granules
- 1 x 1 litre diluter
- 10 x pairs of nitrile gloves
- 10 x aprons
- 10 x clinical waste bags
- 10 x New cardboard scoop and scraper
- 10 x New SoChlor Applicator
- Instructions / Poster
- Recording card
- Wall mount

**SPILLS COVERED INCLUDE:**

- |                        |                           |
|------------------------|---------------------------|
| Blood                  | Cerebrospinal fluid (CSF) |
| Amniotic fluid         | Semen                     |
| Unfixed tissues/organs | Breast Milk               |
| Peritoneal fluid       | Vaginal secretions        |
| Pleural fluid          | Saliva                    |
| Pericardial fluid      |                           |
| Synovial fluid         |                           |



## Biohazard Spill Kit Max

**Code: MJZ555** **All Biohazard Spills** **LP Sorbent Powder**

For use on up to 15 Biohazard Spills (All spills covered).

**NEW**

**KIT CONTENTS:**

- 15 x pair of nitrile gloves
- 15 x apron
- 15 x New LP Sorbent Powder Sachet
- 1 x smart diluter
- 15 x New cardboard scoop and scraper
- 15 x clinical waste bag
- 15 x SoChlor TAB 0.5g x4 sachet packs
- 1 x New SoChlor Applicator
- Instructions / Poster
- Recording card
- Wall mount

**SPILLS COVERED INCLUDE:**

- |                           |                    |
|---------------------------|--------------------|
| <b>Blood</b>              | Semen              |
| Amniotic fluid            | Breast Milk        |
| Unfixed tissues/organs    | Vaginal secretions |
| Peritoneal fluid          | Saliva             |
| Pleural fluid             | <b>Urine</b>       |
| Pericardial fluid         | <b>Vomit</b>       |
| Synovial fluid            | Faeces             |
| Cerebrospinal fluid (CSF) | Sputum             |



Available 2025

## Watch Video Guides

Illustrated and easy to follow video guides for all our Spill Kits & Packs.

Visit: [gvhealth.com/knowledge-centre](http://gvhealth.com/knowledge-centre) or simply scan the QR codes



Biohazard Spill Pack (MJZ003)



UV Spill Pack (MJZ004)



Midi Biohazard Spill Kit (MJZ018)

Biohazard Spill Kit Mini: 2 Pack (MJZ016)

# Cytotoxic Spill Pack

**Code: MJZ015** **LP Sorbent Powder**

Cytotoxic spill packs provide a rapid, safe way to deal with spillages of cytotoxic drugs. Used mainly in the treatment of cancers, cytotoxic drugs are handled in outpatient areas, wards and pharmacies and regularly carried between these locations.

- PACK CONTENTS**
- 1 x New Sorbent Powder Sachet
  - 1 x New cardboard scoop and scraper
  - 1 x anti-splash face mask
  - 1 x apron\*
  - 1 x pair of nitrile gloves\*
  - 1 x cytotoxic waste bag
  - 1 x 10g Cytali Powder Sachet
  - 1 x New Applicator
  - 1 x 250ml diluter

\*ASTM D6978-05: Standard Practice for Assessment of Resistance of Medical Gloves to Permeation by Chemotherapy Drugs



Single use, liquid spill, on hard surface

# Cytotoxic Spill Kit

**Code: MJZ000** **LP Sorbent Powder**

Lightweight and portable for multiple and/or large scale spillages. Our new LP Sorbent Powder absorbs polar non-aqueous spills such as cytotoxic liquids. Absorbent pads and sterile water mop up powder spills. Our custom Cytali solution, designed to remove cytotoxic residue of spill site and surrounding area. Warning signs and colour coded, sealable bags for safe disposal.

**KIT CONTENTS:**

- |                           |                                     |                         |
|---------------------------|-------------------------------------|-------------------------|
| <b>Personal Equipment</b> | <b>Liquid Clean up</b>              | <b>Waste disposal</b>   |
| 3 x nitrile gloves*       | 3 x New Sorbent Powder Sachet       | 1 x instruction poster  |
| 3 x gown / apron          | 3 x New cardboard scoop and scraper | 3 x Cytotoxic waste bag |
| 3 x FFP3 mask             | 1 x tweezers                        | 3 x bag tie             |
| 3 x eye visor             |                                     | 1 x warning sign        |
| 3 x pairs of over sleeves |                                     | 1 x record card         |
| 3 x pairs overshoes       |                                     |                         |
| 1 x kneeling pad          |                                     |                         |
|                           | <b>Decontamination</b>              |                         |
|                           | 3 x 10g Cytali sachet               |                         |
|                           | 1 x 250ml diluter                   |                         |
| <b>Powder Clean up</b>    | 3 x New Applicator                  |                         |
| 3 x sterile water pods    |                                     |                         |
| 3 x absorbent pads        |                                     |                         |



Multi-use, liquid & powder spills, on hard & soft surfaces

\*ASTM D6978-05: Standard Practice for Assessment of Resistance of Medical Gloves to Permeation by Chemotherapy Drugs

# Cytotoxic Waste Bags

UN approved and ADR compliant (18 01 08\*) – our purple cytotoxic waste bags are used for the disposal of infectious waste that is contaminated with cytotoxic or cytostatic materials including chemotherapeutics.

Manufactured with all our customers in mind, infectious waste bags are available in a wide range of specifications, all passing UN packaging tests, ensuring optimum strength and quality.

As well as taking steps to improve sustainability, our Polycare range is a cost-effective and innovative solution to high-quality products, made from advanced polythene and blend technology, so you can be confident that your environment is protected.

- Code: MVN453**
- Product Specifications:**
- Dimensions:** 180 x 711 x 965mm (7 x 28 x 38 inch)
- Total Qty:** Box of 250 (5 Rolls of 50)



# Chemical Spill Pack

**Code: MJZ019**

A single use spill pack for protection against hazardous chemicals and contaminating fluids. The pack absorbs up to 1 litre of chemical spill and is safe and efficient for cleaning up chemical spillages in a laboratory environment. The small pack makes it easy to keep ready for use in all areas and ensure accessibility where it may be needed. The one pack use format prevents fragmentation of components and ensures all items are present when opened.

- PACK CONTENTS:**
- 1 x pair of nitrile gloves
  - 1 x apron / thumb loop gown
  - 1 x face mask/eye shield
  - 2 x surface wipes
  - 1 x disposal bag
  - 2 x absorbent pads
  - 1 x set of tweezers
- Easy to follow illustrated instructions



Scan to watch chemical spill procedure video



# Chemical Spill Kit - 3 pack

**Code: MJZK19** **Refill Pack Format**



Multi-use Spill Kit for fast containment and removal of hazardous chemical spills. The Chemical Spill Kit contains 3 single use packs; each pack absorbs up to 1 litre of chemical spill. The absorbent pad system is tested and approved to be used with over 60 different chemical spills. Each pack also contains all PPE to ensure fast user safety in the event of a hazardous chemical spill. The kit also include a wall bracket to ensure accessibility in areas where it may be needed. The one pack use format prevents fragmentation of components and ensures all items are present when opened.

- Hard Case, Wall Bracket & x3 MJZ019 Packs**
- MJZ019 Pack contents:**
- 1 x pair of nitrile gloves
  - 1 x apron / thumb loop gown
  - 1 x face mask/eye shield
  - 2 x surface wipes
  - 1 x disposal bag
  - 2 x absorbent pads
  - 1 x set of tweezers
- Easy to follow illustrated instructions



Scan to watch chemical spill procedure video



**MJZ019 & MJZK19 fully tested for use on the following chemicals:**

- |                       |                         |                           |                           |
|-----------------------|-------------------------|---------------------------|---------------------------|
| Acetic acid glacial   | Cupric chloride         | Jet fuel (JP-5)           | Sodium hydroxide (50%) ** |
| Acetic acid           | Cyclohexanone           | Kerosene                  | Sodium hypochlorite       |
| Acetone               | Dichloromethane         | Methanol                  | Styrene                   |
| Acetonitrile          | Diethylamine            | Methyl ethyl ketone       | Sulfuric acid (50%)       |
| Aluminum salts        | Dimethylformamide       | Methyl isobutyl ketone    | Sulfuric acid (98%)*      |
| Ammonium fluoride     | Ethanol                 | Mineral oil               | Tetrachloroethylene       |
| Ammonium hydroxide    | Ethyl acetate           | Mineral spirits           | Tetrahydrofuran           |
| Aqueous ammonia (29%) | Formaldehyde            | Naphtha                   | Thionyl chloride          |
| Barium salts          | Gasoline                | Nitric acid (70%)         | Toluene                   |
| Benzyl alcohol        | Glycol ether            | Nitric acid (fuming, 90%) | 1, 1, 1-Trichloroethane   |
| Boric acid            | Hexane                  | Nitrobenzene              | Trichloroethylene         |
| Butanol               | Hydrochloric acid (37%) | Perchloroethylene         | Triethylamine water       |
| Butyl acetate         | Hydrogen peroxide (30%) | Phenol                    | Turpentine                |
| Calcium chloride      | Hydrogen peroxide (50%) | Phosphoric acid (86.7%)   | Xylene                    |
| Carbon disulfide      | Hydrofluoric acid (48%) | Propylene glycol          |                           |
| Carbon tetrachloride  | Isopentyl acetate       | Sodium hydroxide (30%)    |                           |
| Chloroform            | Isopropanol             | Sodium hydroxide (40%)*   |                           |

## High Absorbency Granules

**Code: GVAG200** **UV Spills & Containment**

High Absorbency Granules are designed for spill containment and removal of biohazards such as urine and vomit. In addition to spill containment and removal, granules can be used to ensure safe transportation of liquid clinical waste. High Absorbency Granules are ideal for use with bedpans, vomit bowls and urine bottles.

### APPLICATIONS

- Absorbs body fluids and converts them to a firm non-spill gel.
- An absorbent powder for the prevention of spillages.
- Ideal for use with bedpans, vomit bowls and urine bottles.
- Assists safe transportation of liquid clinical waste.

### SPILLS COVERED INCLUDE:

- Urine
- Vomit
- Faeces
- Saliva
- Sputum



## UV Spill Kit Multi

**Code: MJZ020**

The Urine and Vomit Spill Kit is designed for the safe clean up and disinfection of 6 spillages of urine or vomit. The kit provides all required PPE and specialist products including High Absorbency Granules, IPA disinfectant wipes and GV Neutralising Spray to ensure best practice removal and disinfection of a spill. GV Health's Neutralising Spray removes residual odour where urine or vomit spills have occurred. The spray can be used on hard and soft surfaces.

### KIT CONTENTS

- 1 x 200g High Absorbency Granules
- 1 x 250ml GV Neutralising Spray
- 6 x pairs of nitrile gloves
- 6 x packs of disinfectant surface wipes
- 6 x New cardboard scoop and scraper
- 6 x clinical waste bags
- 6 x aprons
- 1 x instruction poster
- 1 x wall bracket
- 1 x blue hard case

### SPILLS COVERED INCLUDE:

- Urine
- Vomit
- Faeces
- Saliva
- Sputum



## UV Spill Kit Multi Pack Refill

**Code: MJZ021**

### KIT CONTENTS

- 1 x 200g High Absorbency Granules
- 1 x 250ml GV Neutralising Spray
- 6 x pairs of nitrile gloves
- 6 x packs of disinfectant surface wipes
- 6 x New cardboard scoop and scraper
- 6 x clinical waste bags
- 6 x aprons
- 1 x instruction poster

### SPILLS COVERED INCLUDE:

- Urine
- Vomit
- Faeces
- Saliva
- Sputum



## SoChlor Granules

**Code: MFB253**

Designed for use on biohazard / blood or blood stains, body fluids and other secretions spillages. SoChlor 50% NaDCC Granules, absorb, contain and reduce risk of cross contamination or spread of infection

### SPILLS COVERED INCLUDE:

- |                        |                     |
|------------------------|---------------------|
| Blood                  | Cerebrospinal fluid |
| Amniotic fluid         | (CSF)               |
| Unfixed tissues/organs | Semen               |
| Peritoneal fluid       | Breast Milk         |
| Pleural fluid          | Vaginal secretions  |
| Pericardial fluid      | Saliva              |
| Synovial fluid         |                     |



Scan to watch spill kit procedure videos



## Life Protected Sorbent Powder

**Code: LPS350** **All Biohazard Spills** **Cytotoxic Spills**

With superior containment properties and a hazard-free formula, LP Sorbent Granules enable safe usage on a range of spills and surfaces. See full details - page 32

Available separately 2025



# PATIENT INTERACTION



### Aprons and Tabards

Our Aprons and Tabards ensure the highest standards for patient safety.

Manufactured for high performance and developed to be more sustainable.

UKCA marked under UK MDR (2002) and registered with the MHRA. Self-certified as non-sterile Medical Devices (Class I).

# HEALTHCARE APRONS

Our single use Healthcare Aprons are self-certified as non-sterile Medical Devices (Class I), UKCA marked under UK MDR (2002) and registered with the MHRA

## Disposable Aprons

### Overview

Manufactured from high quality material, our range of latex-free aprons provide advanced barrier protection from direct contact with solids or liquids.

Accessible in a wide variety of sizes, colours and lengths – the versatility of these adjustable aprons make them suitable for use across different sectors including healthcare, food preparation and cleaning.

With sustainability at our forefront, selected products from our Polycare range are crafted with a tailored blend of virgin and recycled plastic, PE and PP. Our manufacturing process optimises the use of polymers, replacing the use of virgin polymers with recycled polymers where possible, whilst maintaining durability.

All colours available:



OUR APRONS MANUFACTURED FROM 2024 WILL CONTAIN 30% OR MORE RECYCLED CONTENT



**WHITE**  
Code: **AWR4616 / APW272**  
27 x 46 inch, 16 Micron, Roll of 200



**GREEN**  
Code: **APG109**  
27 x 46 inch, 16 Micron, Roll of 200



**GREEN MD**  
Code: **APG078**  
27 x 53 inch, 35 Micron, Roll of 75



**GREEN HD**  
Code: **APG069**  
34 x 59 inch, 47 Micron, Pack of 25



**RED**  
Code: **APR108**  
27 x 46 inch, 16 Micron, Roll of 200



**CLEAR**  
Code: **ACR4616**  
27 x 46 inch, 16 Micron, Roll of 200



**WHITE HD APRON**  
Code: **APW079**  
27 x 53 inch, 35 Micron, Roll of 75



**PINK**  
Code: **APP112**  
27 x 46 inch, 16 Micron, Roll of 200



**YELLOW**  
Code: **APY110**  
27 x 46 inch, 16 Micron, Roll of 200



**BLUE**  
Code: **APB114**  
27 x 46 inch, 16 Micron, Roll of 200

## Paediatric Apron

### Overview

Our paediatric aprons are an upgraded design from standard polymer aprons. Tailored for specialist use in paediatric wards, their double-layer design provides an optimal barrier for patient interactions and maintains a hygienic environment for paediatric care.

- White paper layer – crafted with absorbent material for the absorption of fluids and spills.
- Plastic laminated backing – acting as a barrier to avoid soiling of garments.



**WHITE PAEDIATRIC APRON**  
Code: **BTB029**  
27.5 x 56 inch, Pack of 250







**CASE STUDY**  
**Patient Safety**

Following preliminary investigations, a 3 month trial at Colchester General Hospital using red 'do not disturb' tabards was undertaken.\* The main findings of this trial were:

- Using tabards resulted in a 71% reduction in interruptions.
- Reducing interruptions lead to a significant decrease in the number of medication incidents.
- The reduction in time taken for drug rounds resulted in patients receiving their medication in a more timely manner.
- As a result of taking less time on drug rounds trained nurses had more time to undertake other nursing activities.

\*Source Nursing Times article: "The effectiveness of drug round tabards in reducing incidence of medication errors."



# DRUG ROUND TABARDS

Single use tabards for nurses on drug rounds. For use in all types of hospital ward including general medical, surgical and SCBU. Also appropriate for use in residential care and nursing homes.

Research by the National Patient Safety Agency also revealed that although about 80% of mistakes caused no harm in 15% of cases there was low harm and in around 5% moderate or severe harm to patients was caused. This equates to a total of over 2,000 moderate or severe cases of harm to patients through drug errors in a year.

GV Health single use tabards were developed in collaboration with a group of NHS Hospital Trusts. This input has directed the overall tabard design and construction as well as colour, wording and fit. Single size tabard adjusts to fit all. Economical – disposable tabard does not require laundering; suitable for short duration tasks.

**Features**

- Can reduce interruptions by over 70%
- Improves Patient Safety by reducing medication administration errors
- Made from 30% recycled polymer blend
- Reduces cost & ensures better hygiene vs laundered Tabards

Order Code	Description	Packing	Units per Roll/Pack	Rolls/Packs per case	Units per case
BT8000	Single use Drug Round Tabards	Pack	125	2	250



Our single use Tabards are self-certified as non-sterile Medical Devices (Class I), UKCA marked under UK MDR (2002) and registered with the MHRA.





**Polycare**

# SPECIMEN BAGS

Specimen bags are used to safely and hygienically transport tissue, blood, urine and stool samples to pathology laboratories for examination. This includes conveyance into chemical pathology, haematology, histopathology and medical microbiology laboratories.

Keep samples contained with a grip-seal closure, and include respective documentation or reports in the attached document wallet for safe transportation. The transparency of our specimen bags allow for barcode reading through the bag, but are also available in a range of colours and prints.

This product should not be used for anything other than tissue, blood, urine and stool samples to be sent for laboratory analysis. Any documentation attached must be disposed of as confidential waste.

#### Features:

- Options on types of seals – no leakage ensured by pillow seal
- All capacities and sizes available
- Advanced print to aid visual confirmation of waste streams
- Flexibility of HDPE & LDPE, regnan and virgin blends used to produce the best product
- Reduced packaging due to lean manufacturing
- Tint colouring to enable bar code reading through the bag
- Mainly made from clear polythene for clear identity but a variety of colours and print available
- Keep samples and reports securely together using Grip seal specimen bags with attached separate document wallet and grip seal for safe transportation



Order Code	Product description	Open width (mm)	Length (mm)	Seal Type	Unit of Issue
<b>SBWP215</b>	Specimen Bags - Clear Resealable with Pouch - Printed in Red	150	230	Pillow	100 per pack
<b>SBWP225</b>	Specimen Bags - Clear Resealable with Pouch - Unprinted	150	230	Pillow	100 per pack
<b>SBWP205</b>	Specimen Bags - Clear Resealable with Pouch - Unprinted	150	140	Pillow	100 per pack
<b>SBR14</b>	Specimen Bag - Clear Resealable with Top Skirt - Unprinted	254	356	Pillow	100 per pack
<b>SBR11</b>	Specimen Bag - Clear Resealable with White Writing Panel	203	280	Pillow	100 per pack
<b>SBWPBIO205</b>	Specimen Bag - Clear Resealable with Pouch Printed with Yellow Background and Red Biohazard Sign	152	140	Pillow	100 per pack
<b>SBWP200</b>	Specimen Bag - Clear Resealable with Pouch - Unprinted	140	140	Pillow	100 per pack

# HEALTHCARE WASTE BAGS



## Next Generation Waste Bags for Healthcare & Public Health

Making the best waste bags requires expertise, experience, and an unwavering commitment to quality. Starting from a deep understanding of users' needs we leverage our knowledge of polythene manufacturing and global supply chains to deliver the bags you need at price points you can afford.

Our excellent manufacturing processes allow us to consistently provide products that excel in quality so that wherever they are used, you can have confidence in their performance.



UN 3291  
CLINICAL WASTE  
unspecified waste

7kg  
MAXIMUM WEIGHT  
Do not fill over 3/4 full.

UN 5H4/Y12/S/20/GB/6968



SCAN TO LEARN MORE ABOUT WASTE STREAMS

This bag has NOT PASSED the tests for bulk transport defined in 7.3.2.6.1(c) / (d) of 4.3.2.4-2.3 (IMDG Code). Not to be used for the transport of Radioactive Material.

# ORANGE CLINICAL BAGS

## Alternative Treatment or Incineration

More sustainable while still ensuring trusted high performance

UN approved and ADR compliant (18 01 03\*) – our orange clinical waste bags are used for the disposal of infectious waste that does not contain chemicals or pharmaceutical medicines.

Manufactured with all our customers in mind, infectious waste bags are available in a wide range of specifications, all passing UN packaging tests, ensuring optimum strength and quality. As well as taking steps to improve sustainability, our Polycare range is a cost-effective and innovative solution to high-quality products, made from advanced polythene and blend technology, so you can be confident that your environment is protected.



Order Code	Waste Stream	Duty UN Max Approved Weight	Size	Dimensions	Packing Details	Bags per Box
<b>MEDIUM DUTY</b>						
<b>OLD5108</b>	Orange Clinical	Medium Duty UN approved weight 5kg	Small	430 x 635 mm 17 x 25 inches	50 per roll 20 rolls in box	1000
<b>ORCP215/800</b>	Orange Clinical	Medium Duty UN approved weight 5kg	Small	285 x 560 x 635 mm 11 x 22 x 25 inches	50 per roll 16 rolls in box	800
<b>ORCP5/500</b>	Orange Clinical	Medium Duty UN approved weight 5kg	Small	220 x 430 x 635 mm 9 x 17 x 25 inches	25 per roll 20 rolls in box	500
<b>ORCP6</b>	Orange Clinical	Medium Duty UN approved weight 6kg	Small	285 x 560 x 635 mm 11 x 22 x 25 inches	25 per roll 16 rolls in box	400
<b>ORCP285</b>	Orange Clinical	Medium Duty UN approved weight 5kg	Large	711 x 965 mm 28 x 38 inches	25 per roll 18 rolls in box	450
<b>ORCP288</b>	Orange Clinical	Medium Duty UN approved weight 8kg	Large	711 x 980 mm 28 x 39 inches	50 per roll 10 rolls in box	500
<b>ORCP8/500</b>	Orange Clinical	Medium Duty UN approved weight 8kg	Large	360 x 711 x 965 mm 14 x 28 x 38 inches	25 per roll 20 rolls in box	500
<b>ORCP258</b>	Orange Clinical	Medium Duty UN approved weight 8kg	Large	325 x 635 x 965 mm 13 x 25 x 38 inches	25 per roll 15 rolls in box	375
<b>ORCP368</b>	Orange Clinical	Medium Duty UN approved weight 8kg	Large	325 x 635 x 915 mm 13 x 25 x 36 inches	25 per roll 8 rolls in box	200
<b>HEAVY DUTY</b>						
<b>ORCP10</b>	Orange Clinical	Heavy Duty UN approved weight 10kg	Large	360 x 711 x 965 mm 14 x 28 x 38 inches	25 per roll 10 rolls in box	250
<b>ORCP2812</b>	Orange Clinical	Heavy Duty UN approved weight 12kg	Large	711 x 980 mm 29 x 38 inches	50 per roll 10 rolls in box	500
<b>BULK ADR 7.3.2.6.2</b>						
<b>ORCP2811</b>	Orange Clinical	Heavy Duty Bulk UN approved weight 11kg	Large	711 x 965 mm 28 x 38 inches	25 per roll 10 rolls in box	250
<b>ORCP2825</b>	Orange Clinical	Heavy Duty Bulk UN approved weight 25kg	Large	711 x 990 mm 28 x 39 inches	10 per roll 10 rolls in box	100
<b>MATTRESS DISPOSAL ADR 7.3.2.6.2</b>						
<b>ORCP4614</b>	Orange Clinical	Heavy Duty UN approved weight 14kg	Large	940 x 1140 x 2500 mm 37 x 45 x 98 inches	5 per roll 1 roll in box	5

 OUR ORANGE CLINICAL BAGS ARE NOW MADE WITH 30% OR MORE RECYCLED CONTENT



**Recycled content**

With sustainability at our forefront, selected products from our Polycare range are crafted with a tailored blend of virgin and recycled plastic, PE and PP. Our manufacturing process optimises the use of polymers, replacing the use of virgin polymers with recycled polymers where possible, whilst maintaining durability.



# YELLOW CLINICAL BAGS

## For Incineration

### More sustainable while still ensuring trusted high performance

UN approved and ADR compliant (18 01 03\*) – our yellow clinical waste bags are used for the disposal of infectious waste that is contaminated with chemicals or pharmaceutical medicines.

Manufactured with all our customers in mind, infectious waste bags are available in a wide range of specifications, all passing UN packaging tests, ensuring optimum strength and quality.

As well as taking steps to improve sustainability, our Polycare range is a cost-effective and innovative solution to high-quality products, made from advanced polythene and blend technology, so you can be confident that your environment is protected.



Order Code	Waste Stream	Duty UN Max Approved Weight	Size	Dimensions	Packing Details	Bags per Box
<b>MEDIUM DUTY</b>						
YRCP215	Yellow Clinical	Medium Duty UN approved weight 5kg	Small	285 x 560 x 635 mm 11 x 22 x 25 inches	50 bags per roll 16 rolls in box	800
YLD5108	Yellow Clinical	Medium Duty UN approved weight 5kg	Small	430 x 635 mm 17 x 25 inches	50 bags per roll 20 rolls in box	1000
YRCP258	Yellow Clinical	Medium Duty UN approved weight 8kg	Large	325 x 635 x 965 mm 13 x 25 x 38 inches	25 bags per roll 15 rolls in box	375
YRCP8	Yellow Clinical	Medium Duty UN approved weight 8kg	Large	360 x 711 x 965 mm 14 x 28 x 38 inches	25 bags per roll 12 rolls in box	300
<b>HEAVY DUTY</b>						
YRCP5	Yellow Clinical	Heavy Duty UN approved weight 5kg	Small	220 x 430 x 635 mm 9 x 17 x 25 inches	25 bags per roll 20 rolls in box	500
YMD10105	Yellow Clinical	Heavy Duty UN approved weight 10kg	Large	711 x 965 mm 28 x 38 inches	25 bags per roll 10 rolls in box	250
YRCP2812	Yellow Clinical	Heavy Duty UN approved weight 12kg	Large	711 x 991 mm 28 x 39 inches	50 bags per roll 10 rolls in box	500
<b>BULK ADR 7.3.2.6.2</b>						
YRCP2811	Yellow Clinical	Heavy Duty Bulk UN approved weight 11kg	Large	711 x 965 mm 28 x 38 inches	25 bags per roll 10 rolls in box	250
YRCP2825	Yellow Clinical	Heavy Duty Bulk UN approved weight 25kg	Large	711 x 991 mm 28 x 39 inches	10 bags per roll 10 in box	100
<b>MATTRESS DISPOSAL ADR 7.3.2.6.2</b>						
YRCP4614	Yellow Clinical	Heavy Duty UN approved weight 14kg	Large	940 x 1140 x 2500 mm 38 x 45 x 99 inches	5 bags per pack 1 pack in box	5



OUR YELLOW CLINICAL BAGS ARE NOW MADE WITH 30% OR MORE RECYCLED CONTENT



### UN Approved

Due to the nature of its contents and its associated risks, clinical waste is classified as dangerous goods. This means there is high regulation around how the waste is stored, transported and disposed of. This primarily applies to infectious chemically-contaminated clinical waste and cytotoxic waste. Bags storing this waste must be tested and approved in accordance with UN regulation to maintain the highest and safest quality.

Clinical waste bags must be able to safely withstand a certain amount of weight for safe carriage and transportation. Our ADR compliant bags have a corresponding UN approved weight limit to indicate maximum approved capacity of the product. This is obtained by putting forward the product for a UN packaging test. This testing assesses the integrity of the product, and if compliant, is approved.

On clinical waste bags that require this approval, the UN number will be printed across the body of the bag to easily identify ADR compliance.



# TIGER STRIPE BAGS

## Offensive Waste

### More sustainable while still ensuring trusted high performance

Our tiger bags are used for the disposal of offensive waste. As the first to manufacture via the coextrusion of two polymers, the yellow and black striping on our offensive waste bags are visible from any angle, avoiding confusion with other yellow clinical waste bags.

Manufactured with all our customers in mind, infectious waste bags are available in a wide range of specifications, all passing UN packaging tests, ensuring optimum strength and quality.

As well as taking steps to improve sustainability, our Polycare range is a cost-effective and innovative solution to high-quality products, made from advanced polythene and blend technology, so you can be confident that your environment is protected.

### Recycled content

With sustainability at our forefront, selected products from our Polycare range are crafted with a tailored blend of virgin and recycled plastic, PE and PP. Our manufacturing process optimises the use of polymers, replacing the use of virgin polymers with recycled polymers where possible, whilst maintaining durability.

Bags made with a blend of polymers create a sustainable approach to waste management, with a robust range of options for different requirements. We aim to find a tailored solution to your needs and budget.



## Offensive Waste

### Can go into offensive waste

- Soiled nappies, sanitary towels and incontinence pads, bed pans
- Non-infectiously contaminated PPE, masks, and dressings (small injuries)
- Empty catheter bags, clean IV bags

### Cannot go into offensive waste

- Liquid offensive waste (liquid faeces, urine) To be disposed of in sluice
- PPE, dressings, masks with infectious contamination
- Chemical, pharmaceutical or cytotoxic contaminated materials



SCAN TO LEARN MORE ABOUT WASTE STREAMS

Order Code	Waste Stream	Duty	Size	Dimensions	Packing Details	Bags per Box
<b>MEDIUM DUTY</b>						
TSLD5101	Offensive Waste	Medium Duty 5kg	Small	432 x 635 mm 17 x 25 inches	50 bags per roll 20 rolls in box	1000
TSMD8110	Offensive Waste	Medium Duty 8kg	Large	711 x 965 mm 28 x 38 inches	25 bags per roll 18 rolls in box	450
<b>HEAVY DUTY</b>						
TSHD1011	Offensive Waste	Heavy Duty 10kg	Large	711 x 990 mm 28 x 39 inches	25 bags per roll 12 rolls in box	300
<b>EXTRA HEAVY DUTY</b>						
GVTS/5	Offensive Waste	Extra Heavy Duty 20kg	Large	457 x 711 x 965 mm 18 x 28 x 38 inches	10 bags per roll 10 rolls in box	100



OUR TIGER BAGS ARE NOW MADE WITH 30% OR MORE RECYCLED CONTENT

# CYTOTOXIC BAGS

## Cytotoxic & Cytostatic Waste

### More sustainable while still ensuring trusted high performance

UN approved and ADR compliant (18 01 08\*) – our purple cytotoxic waste bags are used for the disposal of infectious waste that is contaminated with cytotoxic or cytostatic materials including chemotherapeutics.

Manufactured with all our customers in mind, infectious waste bags are available in a wide range of specifications, all passing UN packaging tests, ensuring optimum strength and quality.

As well as taking steps to improve sustainability, our Polycare range is a cost-effective and innovative solution to high-quality products, made from advanced polythene and blend technology, so you can be confident that your environment is protected.



Order Code	Waste Stream	Duty UN Max Approved Weight	Size	Dimensions	Packing Details	Bags per Box
<b>MEDIUM DUTY</b>						
<b>MVN453</b>	Cytotoxic Clinical	Medium Duty UN approved weight 5kg	Large	711 x 965 mm 28 x 38 inches	50 bags per roll 5 rolls in box	250
<b>PRCP8</b>	Cytotoxic Clinical	Medium Duty UN approved weight 8kg	Large	360 x 711 x 965 mm 14 x 28 x 38 inches	25 bags per roll 12 rolls in box	300

 OUR CYTOTOXIC BAGS ARE NOW MADE WITH 30% OR MORE RECYCLED CONTENT



SCAN TO LEARN MORE ABOUT WASTE STREAMS

## GENERAL WASTE

### More sustainable while still ensuring trusted high performance

Our domestic bags are manufactured for the disposal of municipal waste. Robust and versatile, domestic waste bags can be used in non-clinical and public areas.

As well as taking steps to improve sustainability, our Polycare range is a cost-effective and innovative solution to high-quality products, made from advanced polythene and blend technology, so you can be confident that your environment is protected.



SCAN TO LEARN  
MORE ABOUT  
WASTE STREAMS

## BLACK WASTE BAGS

Order Code	Waste Stream	Product Details	Size	Dimensions	Packing Details	Bags per Box
<b>UNPRINTED</b>						
<b>BRGN4170</b>	Black General	Wheelie Bin Liner	Extra Large	711 x 1042 x 1397 mm 28 x 41 x 55 inches	50 bags per pack 1 pack in box	50
<b>PRINTED WITH "NHS PROPERTY FOR DOMESTIC WASTE ONLY"</b>						
<b>BVDP17</b>	Black General	Black Bag Tested weight 4kg	Small	432 x 635 mm 17 x 25 inches	50 bags per roll 10 rolls in box	500
<b>BVDP21</b>	Black General	Black Bag Tested weight 4kg	Large	535 x 595 mm 21 x 24 inches	50 per roll 10 rolls in box	500
<b>BVDP28L</b>	Black General	Black Bag Tested weight 5kg	Large	711 x 950 mm 28 x 38 inches	50 bags per roll 6 rolls in box	300
<b>BVDP28</b>	Black General	Black Bag Tested weight 5kg	Large	711 x 950 mm 28 x 38 inches	50 bags per roll 6 rolls in box	300

## CLEAR BAGS

Order Code	Waste Stream	Product Details	Size	Dimensions	Packing Details	Bags per Box
<b>UNPRINTED</b>						
<b>CVD295</b>	Clear General	Clear Bag Tested weight 5kg	Large	447 x 737 x 965 mm 18 x 29 x 38 inches	25 bags per roll 20 rolls in box	500
<b>CRN2830</b>	Clear General	Clear Bag Tested weight 5kg	Large	600 x 755 mm 24 x 30 inches	100 bags per roll 5 rolls in box	500
<b>CRN2934</b>	Clear General	Clear Bag Tested weight 5kg	Large	730 x 845 mm 29 x 33 inches	50 bags per roll 5 rolls in box	250
<b>CRGN15241000</b>	Clear General	Square clear bag Tested weight 8kg	Small	610 x 610 mm 24 x 24 inches	50 bags per roll 20 rolls in box	1000
<b>GV035</b>	Clear General	Clear Bag Tested weight 8kg	Large	800 x 1150 mm 32 x 45 inches	50 bags per roll 5 rolls in box	250
<b>GV044</b>	Clear General	Clear Bag Tested weight 10kg	Large	965 x 1067 mm 38 x 42 inches	10 bags per roll 20 rolls in box	200
<b>PRINTED WITH "NHS PROPERTY FOR DOMESTIC WASTE ONLY"</b>						
<b>CVDP21</b>	Clear General	Clear Bag Tested weight 4kg	Small	535 x 595 mm 21 x 24 inches	50 bags per roll 10 rolls in box	500
<b>CVDP28</b>	Clear General	Clear Bag Tested weight 5kg	Large	711 x 960 mm 28 x 38 inches	50 bags per roll 9 rolls in box	450
<b>CENTRIFUGE BAG</b>						
<b>MVB007</b>	N/A	Centrifuge Bag	Small	204 x 308 mm 8 x 12 inches	100 bags per pack 10 packs in box	1000
<b>BLOOD BAG</b>						
<b>MVN801</b>	N/A	Clear Bag	Large	762 x 914 mm 30 x 36 inches	50 bags per roll 4 rolls in box	200



OUR CLEAR AND BLACK BAGS ARE MADE WITH 30% OR MORE RECYCLED CONTENT



# TAGS & CABLE TIES


## Waste Tracking & Identification



Waste tracking and identification is now mandatory throughout the NHS, in independent hospitals and clinics in the UK and in many countries world-wide. Clinical waste and many other categories of healthcare waste must be traceable back to the original source.

GV Health offer a range of purpose designed tapes, tags and ties to help meet the requirement for clinical waste tracking, tracing and identification. Items are pre-printed with your specific tracing codes.

Order Code	Description	Size	Print	Units per Box
<b>TAGS</b>				
<b>GVTLB200SN</b>	Bag Tablock Ties - Yellow - Printed Serial	200 x 4.8mm	Printed	10000
<b>GVZIPPY200X</b>	Bag Ties Plain Cable Ties - Yellow	200 x 4.8mm	Unprinted	1000
<b>GVZIPB200X</b>	Bag Ties Plain Cable Ties - Black	200 x 4.8mm	Unprinted	1000
<b>UVTT300-N</b>	Bag Ties - Flag Head - Black	300mm	Unprinted	4000
<b>GVZIPPY300X</b>	Bag Ties Plain Cable Ties - Yellow	300mm	Unprinted	1000

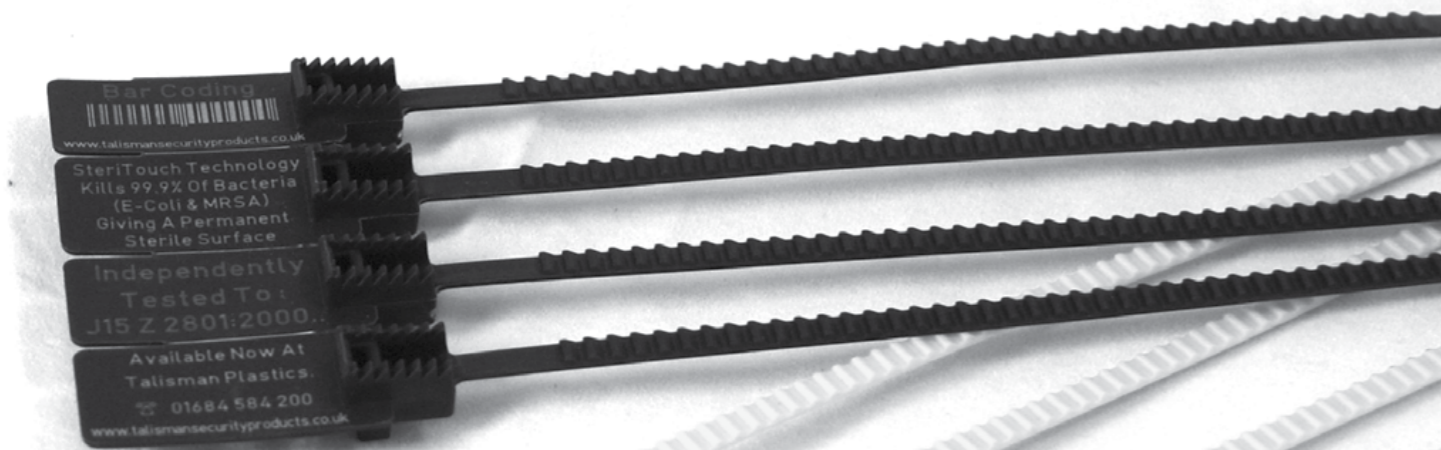

More colours & customisation options available

We can also provide customised items to meet special requirements. Pre-printed cable ties - Identification tags for hospital waste bags. Offer full traceability to the hospital and ward. Flag ties - Supplied in easily separated strips of 10 to aid handling and sequential use. Large flag can be custom printed, allowing full traceability.



## COLOUR WASTE BAGS

### Recycling Sortation



Order Code	Type	Description	Size	Dimensions	Packing Details	Bags per Box
<b>GVGNB</b>	Colour	Green Medium Duty Bag Unprinted 5kg	Large	736 x 940 mm	50 bags per roll	500
<b>GWBL4660</b>	Colour	Green Tint Wheelie Bin Liner Unprinted	Large	876 x 1321 x 1524 mm	50 bags per pack	50

## LAUNDRY BAGS

### Safe handling of infectious laundry

Order Code	Type	Description	Size	Dimensions	Packing Details	Bags per Box
<b>SOLS6950</b>	Laundry	Red Soluble Strip Bag 6kg	Large	695 x 950 mm	50 bags per pack	250

## CONFIDENTIAL WASTE

Order Code	Type	Description	Size	Dimensions	Packing Details	Bags per Box
<b>GV001/GRY</b>	Confidential	Grey printed "Confidential Waste"	Large	483 x 762 x 959 mm	25 bags per roll	250

## NEW

Available for FREE to all users of GV Health products

# WASTE SEGREGATION TRAINING VIDEO



Managing waste generated in Healthcare can seem complex. Our Healthcare Waste Management trailer introduces our informative video guide, covering waste segregation, disposal, cost savings and ADR.

[gvhealth.com/knowledge-centre](http://gvhealth.com/knowledge-centre)



Life. Protected.



ISO 13485  
Certificate No. MD1348515



ISO 9001  
Certificate No. GB2004284



ISO 14001  
Certificate No. EM2001057



GV Health Ltd, Hall House, 2 Arlington Court, Whittle Way,  
Arlington Business Park, Stevenage, Hertfordshire, SG1 2FS, UK  
T +44 (0)1920 463 098 E info@gvhealth.com

[gvhealth.com](http://gvhealth.com)